



*A guide for:*  
**Economy Lift Exam Table**

**SHOR-LINE**  
SCHROER MFG. CO.

**▲WARNING** This product is for animal use only. We do not condone this product for other uses.

**▲WARNING** To reduce the risk of injury, the user must read and understand this guide before installation of this product.

{Save This Guide For Future Reference}

PN: 006.1240.00  
Updated: 2/2006

# Introduction

Thank you for purchasing Shor-Line products. As a leader in animal care equipment, our commitment to you is exactly the same as it was in 1927.

We hope this guide answers any questions you might have in regards to use, care or installation. If you can't find the answers to your questions, have further questions or would like to purchase additional products, please contact us at 800.444.1579 or 913.281.1500. If in Europe phone +44 1446 77 20 41.

Consider keeping this guide in a safe and convenient place for future assistance.

*For your safety, please follow the instructions on these labels:*

## **▲CAUTION ▲WARNING**

---

**▲WARNING** This product is for animal use only. We do not condone this product for other uses.

**▲WARNING** To reduce the risk of injury, the user must read and understand this guide before installation of this product.

---

{Shor-Line}  
Schroer Manufacturing Company  
511 Osage Ave.  
Kansas City, Kansas 66105, USA  
phone 800.444.1579 fax 913.281.5339  
email [guides@shor-line.com](mailto:guides@shor-line.com)  
[www.shor-line.com](http://www.shor-line.com)

{Shor-Line Limited}  
Vale Business Park  
Llandow,  
Vale of Glamorgan CF71 7PF  
United Kingdom  
phone +44 1446 77 20 41 fax +44 1446 77 36 68  
email [quality@shor-line.co.uk](mailto:quality@shor-line.co.uk)  
[www.shor-line.co.uk](http://www.shor-line.co.uk)



## Unpacking & Inspecting Steps

---

### Step 1

Locate the Packing Slip that describes the parts and quantities shipped. The parts for your Economy Lift Exam Table should be carefully examined upon arrival for possible damage in shipment. If any damage is found, notify Shor-Line immediately.

### Step 2

Review Shor-Line's warranty information at: [www.shor-line.com/warranty](http://www.shor-line.com/warranty).

### Step 3

Remove and separate all contents from the shipping containers. Use the Packing Slip to ensure nothing has been lost.

### Step 4

Check content for concealed damage.

*If Concealed Damage:*


Call Shor-Line immediately so we can expedite replacements or repairs.

**IMPORTANT:** After fifteen (15) calendar days of receipt of merchandise, this policy becomes void.

• Economy Lift Exam Table, Electric Lift, 120 V	Item # 903.1240.00
• Economy Lift Exam Table, Electric Lift, 230 V	Item # 903.1240.01
• Economy Lift Exam Table, Hydraulic Lift	Item # 903.1240.02


## Economy Lift Exam Table, Electric Lift

---

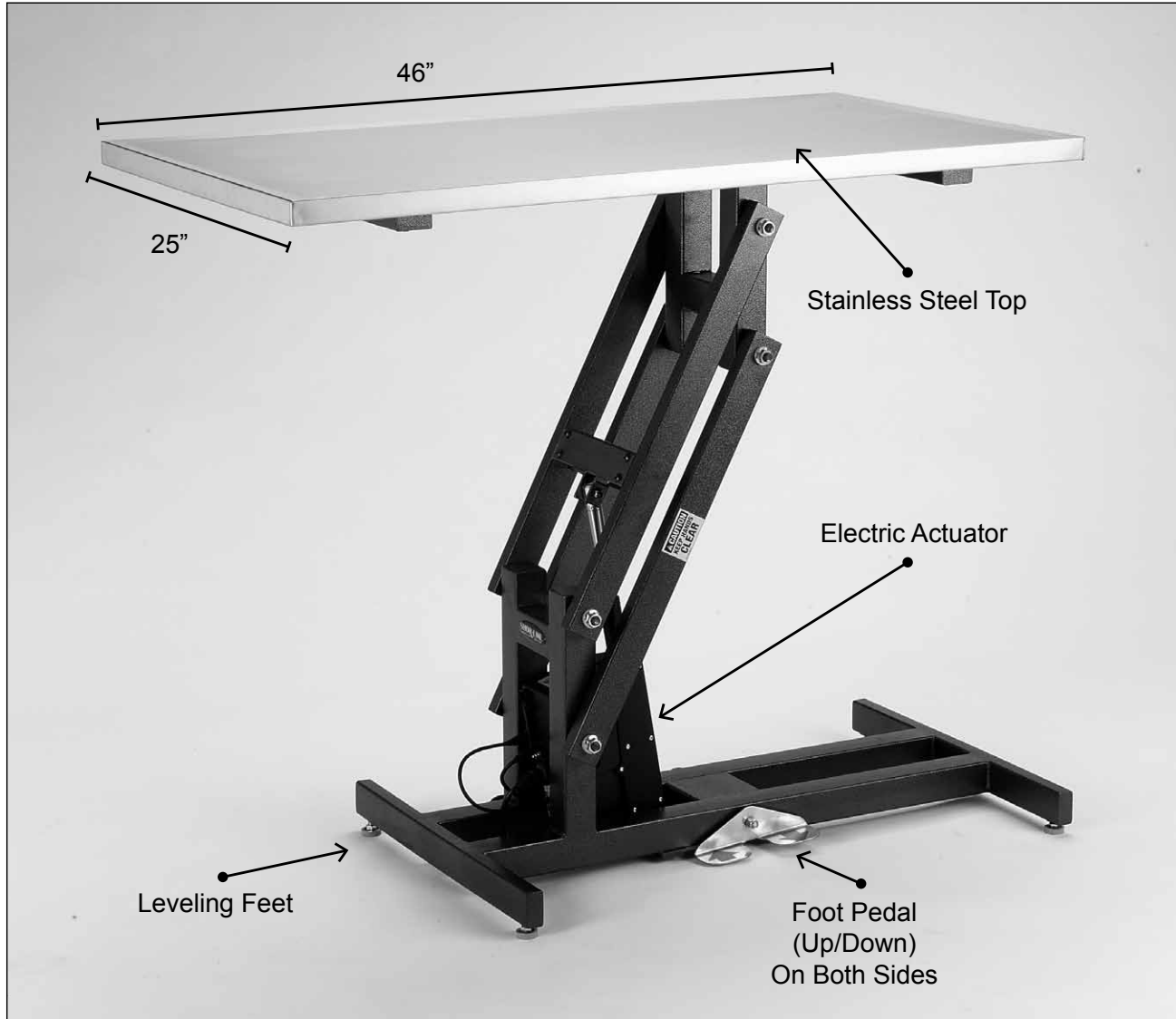

## Economy Lift Exam Table, Hydraulic Lift

---

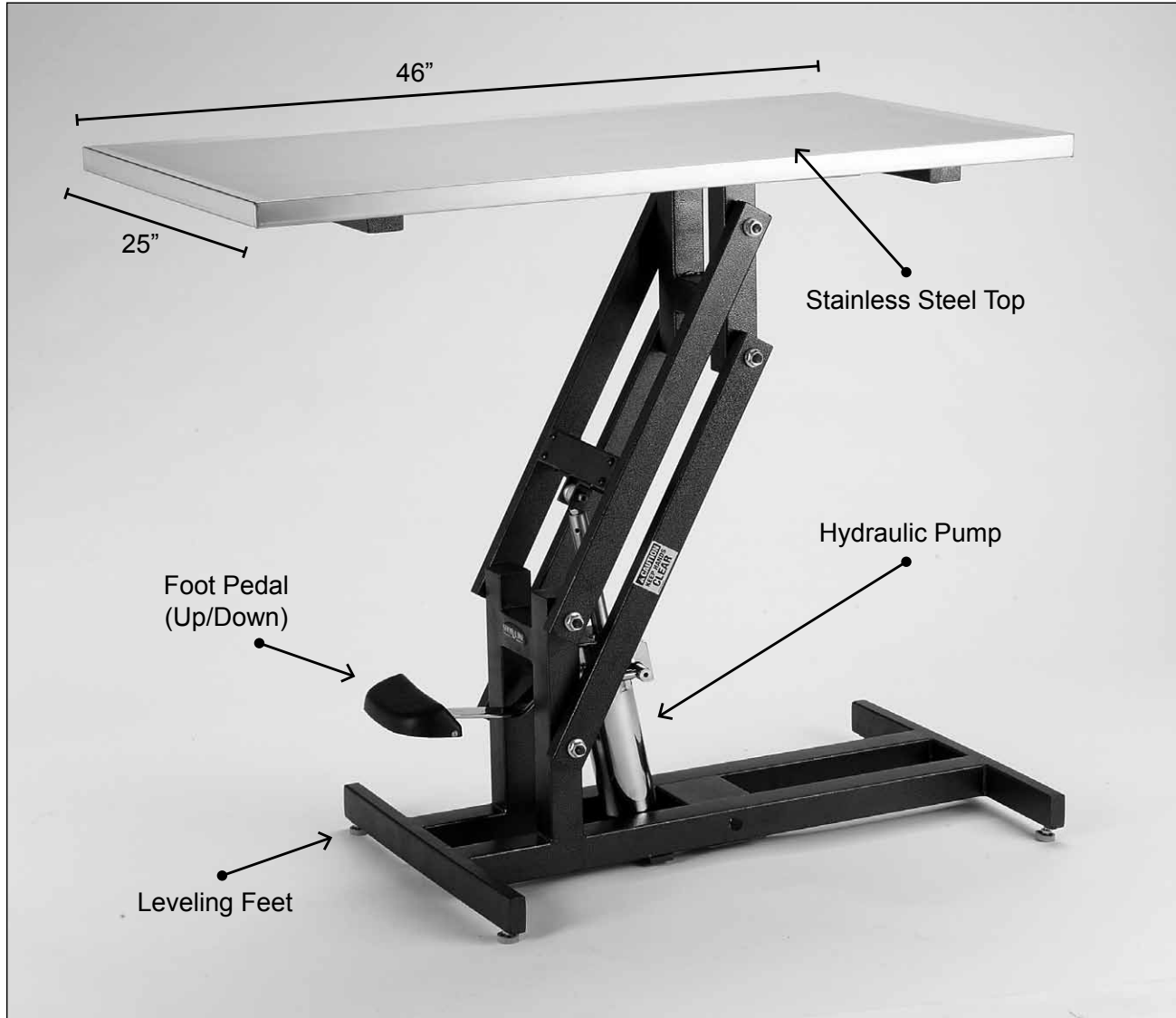

### {To Raise the Table Top}

Push the foot pedal down about 1/2 way from the bottom and release. Continue this until the top is at the desired height.



### {To Lock the Table Top at the Desired Height}

Lift the foot pedal up as far as it will go and release.

### {To Lower the Table Top}

Push the foot pedal down as far as it will go and hold.

Overall Dimensions

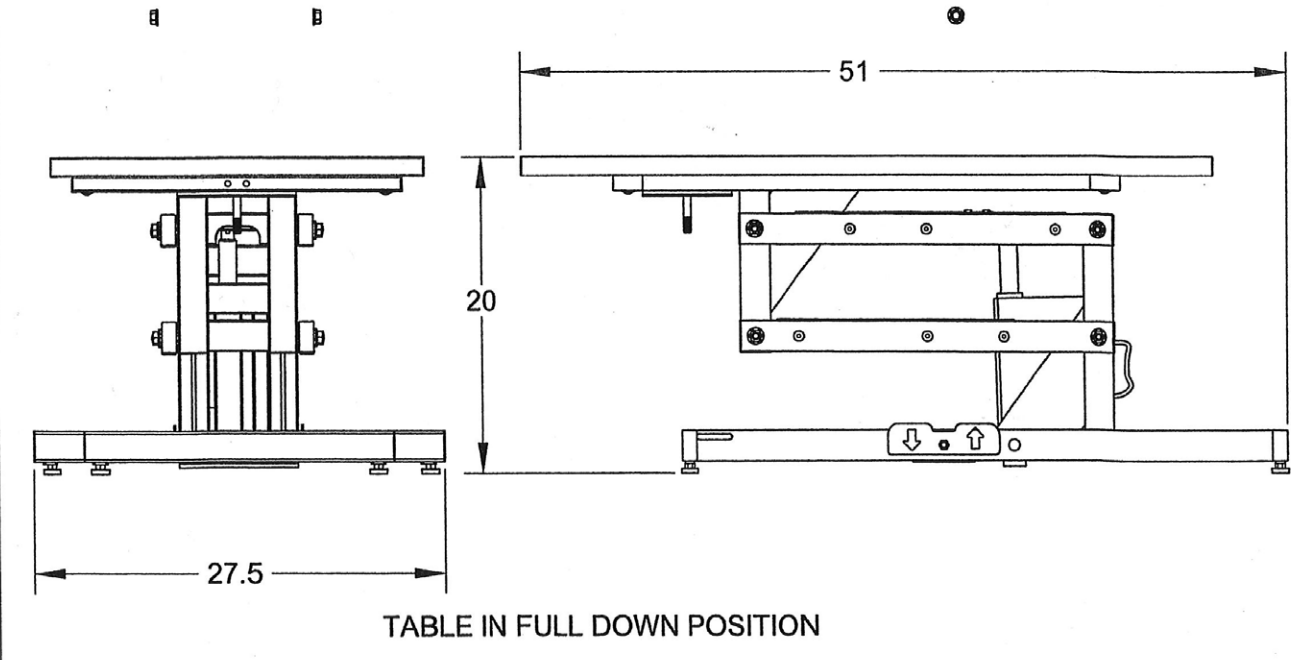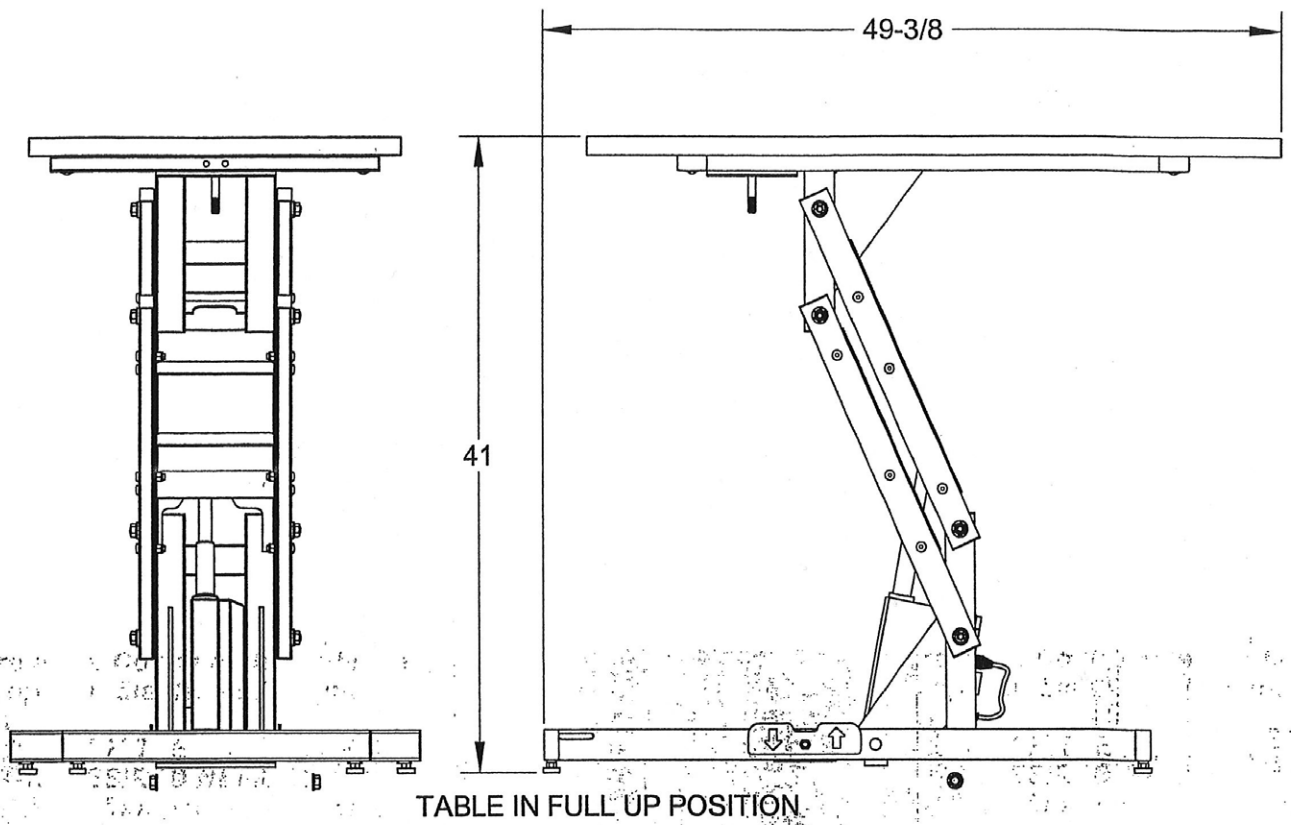

## Specifications

---

### {Economy Lift Exam Table, Electric Lift}


#### *Actuator:*

- 24 Volt DC Motor
- 2.7 AMPS maximum at 250 pound (1112 Newton) load.
- Maximum Shaft Travel - 6 inches (152 mm).
- 120 Volt AC Power Cord for 120 V Model.
- 230 Volt AC Power Cord for 230 V Model.
- DC Control Module Mounted to Actuator.

#### *Table:*

- Maximum Height (Floor to Top of Table) - 41 inches (1041 mm)
- Minimum Height (Floor to Top of Table) - 20 inches (508 mm)
- Time from Minimum to Maximum Height (No Load) - 14 seconds
- Time from Maximum to Minimum Height (No Load) - 13 seconds
- Max Dimensions - Standard Table - 41 inches (1041 mm) High x 51 inches (1295 mm) Long x 27.5 inches (699 mm) wide
- Min Dimensions - Standard Table - 20 inches (508 mm) High x 49.5 inches (1257 mm) Long x 27.5 inches (699 mm) wide

#### **WIRING DIAGRAM-- FOOT SWITCH**

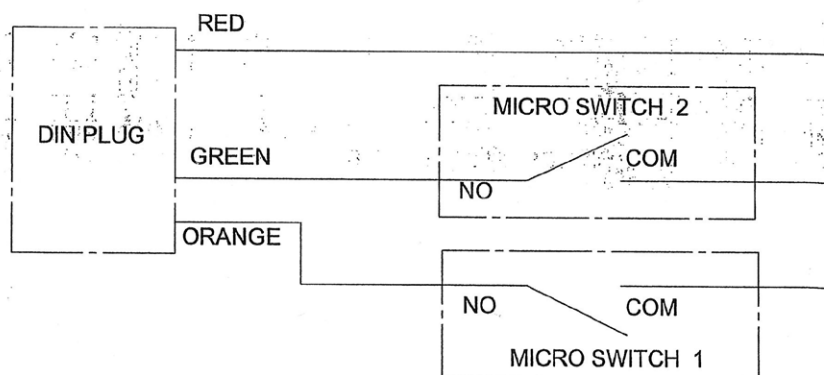

### {Economy Lift Exam Table, Hydraulic Lift}

#### *Actuator:*

- Manual Hydraulic Pump, Foot Pedal Actuated.
- Maximum 250 pound (1112 Newton) Load.
- 

#### *Table:*

- Maximum Height (Floor to Top of Table) - 41 inches (1041 mm)
- Minimum Height (Floor to Top of Table) - 20 inches (508 mm)
- Time from Minimum to Maximum Height (No Load) - 14 seconds
- Time from Maximum to Minimum Height (No Load) - 13 seconds
- Max Dimensions - Standard Table - 41 inches (1041 mm) High x 51 inches (1295 mm) Long x 27.5 inches (699 mm) wide
- Min Dimensions - Standard Table - 20 inches (508 mm) High x 49.5 inches (1257 mm) Long x 27.5 inches (699 mm) wide

# Stainless Steel Care & Cleaning

---

## {Quality & Life of Your Product}

IMPORTANT: Routine maintenance will ensure quality and life of your long-term investment and protect animals from transmittable disease and infection

IMPORTANT: To safeguard your investment and rights to the Shor-Line's Warranty, it is important that you know factors that degrade stainless steel. Corrosion, pitting, and discoloration can occur if:

- Bleach, deodorizing agents, disinfectants, and sanitizers (chlorides, bromides, iodides and thiocyanates) stand wet on your product for a prolonged period (over two minutes).
- Bleach, deodorizing agents, disinfectants, and sanitizers are used frequently without being neutralized with clean water after each application.
- Acidic urine is not removed rapidly and surfaces are not washed thoroughly.
- Animal hair, holding urine and fluids, is allowed to collect around legs, hinges, casters, and welded wire intersections.
- Solutions like liquid bleach and antibacterials are allowed to evaporate and dry on surfaces.
- Particles from steel wool, metal fiber brushes, or metal scrapers (which are never recommended for cleaning stainless) get embedded in your product's surface, rust, and give the false impression that the stainless steel is defective.
- Your water source (local or well) is heavy in iron content and water is allowed to dry on your product.
- High levels of pollution (automotive exhaust, aggressive gases, soot, iron oxides, chemicals, etc.) or chloride (from sea or de-icing salts) exist in the environment and frequent, thorough cleaning (manual or rain-washing) does not occur.

## {Care & Cleaning Best Practices}

IMPORTANT: You play an essential role in safeguarding your investment and rights to Shor-Line's Warranty. We strongly recommend:

- Daily or weekly cleaning of stainless steel to prolong its longevity and value and to help you maintain its finish.
- Using the mildest cleaning procedure that will do the job effectively.
  - Ordinary waste deposits and fluids can usually be removed with soap and water. Rinse thoroughly with clear water and dry completely with a soft cloth to discourage water spotting.
  - Minor scale build up and some hard water spotting can be removed by washing with vinegar followed by a clearwater rinse and thorough drying.
  - If rubbing, brushing, or scraping is required, use only stainless steel wool, brushes with plastic, nylon or stainless steel fibers, and scrubbers made of plastic or nylon. Ordinary steel wool or steel brushes should never be used.
  - Stubborn deposits or tightly adhering debris may require a commercial cleaning product. Scouring powders will scratch stainless and should not be used. Always review label statements on cleaning products regarding use on metals.
  - Cleaning is not completed until fresh water has been used to thoroughly rinse all residue of salt, bleach, deodorizing agents, disinfectants, and sanitizers.
- Rubbing, drying, and polishing should always be in the direction of the product's grain so scratches will blend into the surface. Efforts opposite the grain may produce a scratched surface appearance.
- Rubbing a soft cloth moistened with olive oil, undiluted white or cider vinegar, or club soda may brighten stainless steel and help remove surface streaks and heat stains.
- Never use paint, lacquer, or varnish as a substitute or shortcut for properly caring and cleaning stainless steel.
- Never intentionally stamp or etch stainless surfaces.



## {Protecting Against Disease & Infection}

IMPORTANT: Shor-Line shares your commitment to health and safety. This information is intended help protect your animals, as well as, safeguard your investment.

- Bleaches can be harmful to animals. Bleaches can react with acids and release chlorine gas that is highly toxic.
- Because of its strong corrosive and toxic properties, bleach is diluted with water. The solution's potency (shelf-life) quickly degrades as the active chlorine dissipates with water. Animal care professionals, therefore, cannot know with certainty how effective the bleach solution is over time without implementing a formal monitoring system.
- If chemical deodorizers, disinfectants, or sanitizers are used, protect your animals from respiratory problems and infections by using clean water to wash residue from all surfaces and dry thoroughly.
- Clorox® Company warns: "We cannot recommend Clorox® to clean metal surfaces, including stainless steel. Prolonged or repeated exposure of metal to strong solutions of our product eventually can cause discoloration or corrosion. To reduce the possibility of metal discoloration, we suggest that exposure time to Clorox® be less than five minutes and that surfaces be rinsed well after application of the solution."

# SHOR-LINE®

SCHROER MFG. CO.

toll free	local	fax
800	913	913
444.1579	281.1500	281.5339

Shor-Line/Schroer Manufacturing Co.  
511 Osage Avenue  
Kansas City, KS 66105 U.S.A.

[www.shor-line.com](http://www.shor-line.com)