

OXYGEN THERAPY CONVERSION DOORS

Increased ventilation option for larger cages


Conversion Doors

Benefits:

The polycarbonate Conversion Door pairs with a new-style O2 door (sold separately) to provide an improved, oxygen therapy cage for larger patients. To retrofit, simply remove the grill door and follow kit directions to add the soft, closed-cell foam gasket and placing the door.

Features:

- 3/8"-thick, clear polycarbonate door.
- A soft, closed-cell foam gasket surrounds the door to help retain oxygen levels in the cage.
- Upper and lower adjustable vents allow for reduction in condensation on the interior door surface.
- We placed vents to the side to provide a good view of the pet.
- **NOTE:** If retrofitting from a Shor-Line Type I or Type II cage, you will need an additional retrofit kit.

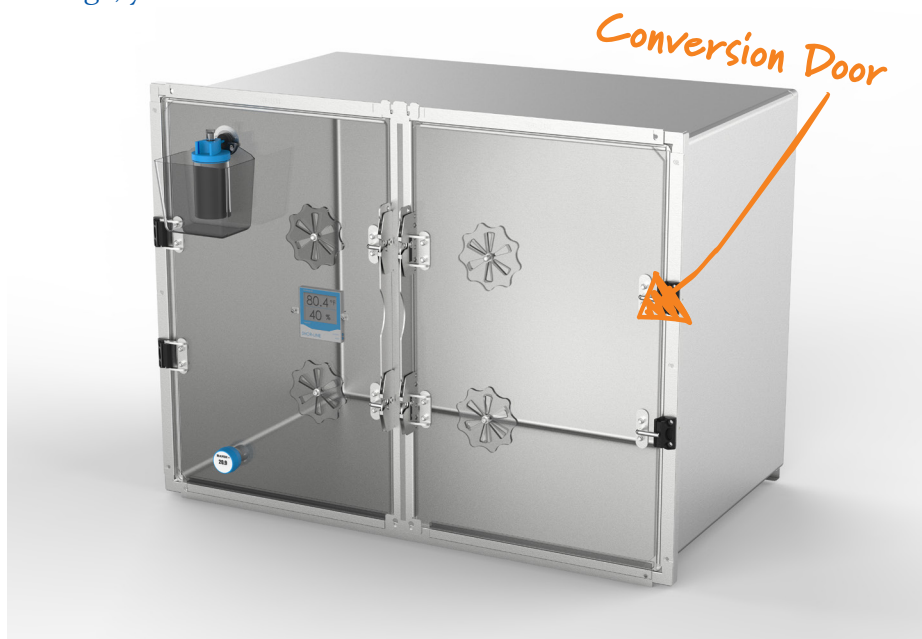

The Conversion Door, on the right, pairs with an Oxygen Therapy Door, sold separately on left, to create a double door O2 cage.

Cage, O2 door and O2 analyzer sold separately

914.2430.21 24"W x 30"H Single Conversion Door 35 lbs.

914.2436.21 24"W x 36"H Single Conversion Door 43 lbs.