

PREFERRED LIFT EXAM TABLE

This may be the best floor-mounted wall flush lift table you will ever use.


Benefits:

This floor-mounted lift table becomes like another pair of hands. The latest technology lets you adjust the height with a key fob or a wireless foot switch. The ergonomic design of the Unified Top allows you to remain face-to-face with pet owners, reducing stress for everyone.

The stylish Stainless Steel Unified Top has a unique backboard design with a tiedown bar that allows multiple restraint positions to handle fractious patients.

Features:

- Patented Stainless Steel exam top has a radiused design that provides space for the pet as well as laptops or other notetaking items. U.S. Patent #D742,526.
- The radiused corners provide additional safety.
- The “T” shape lets you position the pet parallel to the backboard or perpendicular to the backboard with the pet facing outward.
- The spacious and multi-position table top allows the pet owner and veterinarian to face each other for a stress-free and secure exam.
- Available with blue or black edge banding.
- State-of-the-art lift effortlessly and quietly moves the table along a sleek back column using precision linear guides.
- The low-stress actuator lifts the patent from 7¾” to 36”H (counter height), and is virtually maintenance free.
- Helps save your back by both lifting the pet and allowing each veterinarian to easily select the best examination height.
- For solo practitioners or treatment area exams, the backboard and tiedown help you position the animal so both hands are available for examining, securing and soothing the pet. There also are two cam cleats in the back for extra security.
- The unit mounts flush against most any wall. It secures to the floor and can be further mounted to studs in the wall—but not required.
- Each unit comes with two key fobs (battery operated), a rechargeable, wireless foot switch, a foot switch charger and a hard-wired back-up switch on the table itself.

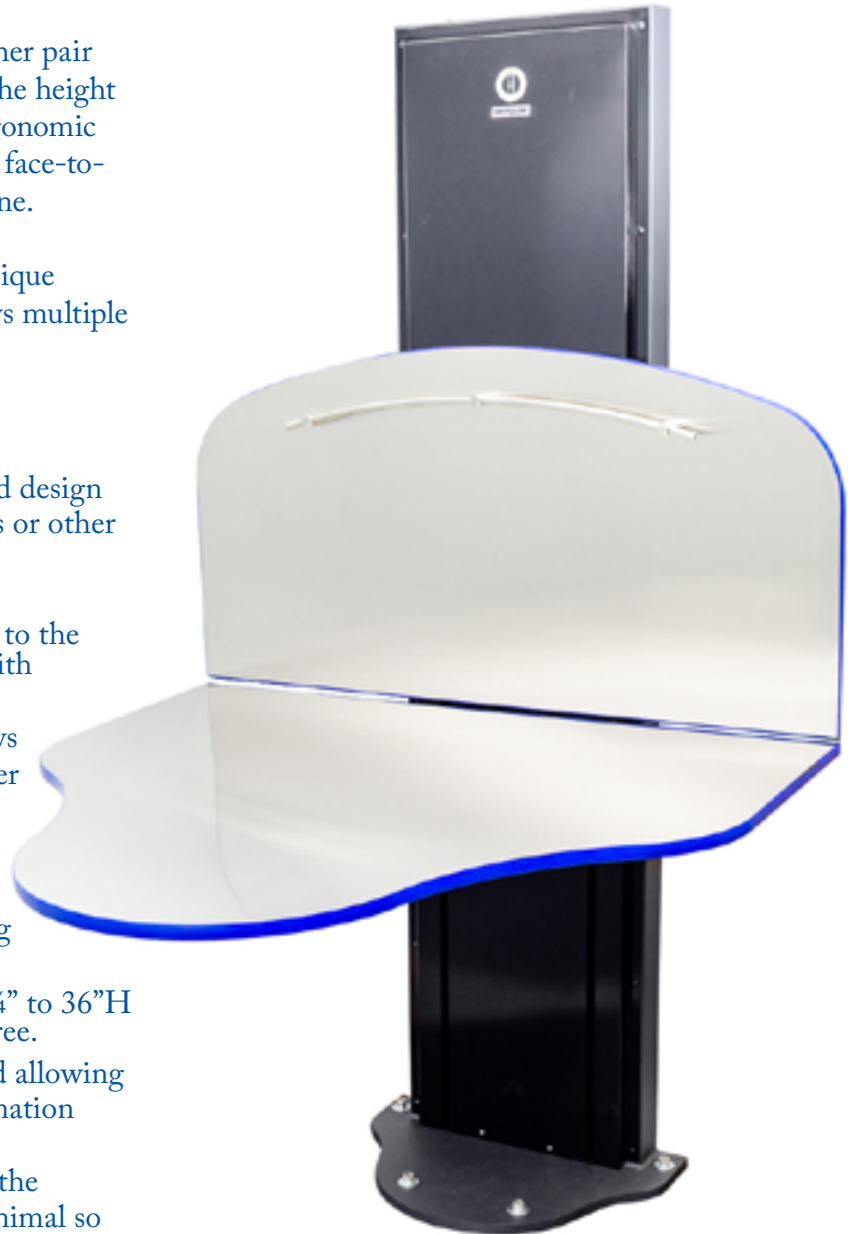

903.3400.00 - Black edge
903.3400.01 - Blue edge

Specs

- Lifts up to 225 lbs.
- Table is 44”W by 41¼”L at its longest point.
- When mounted against the wall, the unit’s column is 75”H and the table extends out 48¼” at its longest point.
- The Stainless Steel backboard is 24-7/8”H.